

Gel Blot Analyzer v9.3

User Manual

Date: September 4, 2026

Anindya Karmaker (Developer)

Email: anindyakarmaker1996@gmail.com and akarmaker@ucdavis.edu

Table of Contents

1. Introduction.....	3
2. Getting Started.....	3
3. Core Concepts.....	10
4. Basic Operations: File & View.....	11
5. Image Adjustments & Transformations.....	15
6. Annotations: Markers & Shapes.....	19
7. Overlap Images.....	23
8. Quantitative Analysis Workflow.....	26
9. Molecular Weight Prediction.....	32
10. Advanced Features.....	34
11. Keyboard Shortcuts.....	41
12. Troubleshooting.....	43

1. Introduction

1.1 What is Gel Blot Analyzer (GBA)?

Gel Blot Analyzer is a desktop application designed for the analysis of images from gel electrophoresis, Western blots, and similar assays. It provides a comprehensive suite of tools for image annotation, lane quantification, distance quantification, protein and glycan analysis with molecular weight prediction in a single, user-friendly interface.

Key Features

- **Advanced Annotations:** Place standard molecular weight markers automatically or manually, custom text labels, arrows, lines, and rectangles with full control over font, size, and color.
- **Gel Analysis:** Adjust the black point, white point and/or gamma to visualize the protein bands clearly.
- **Powerful Quantification:** Define lane areas using precise quadrilaterals (for skewed lanes) or rectangles. Analyze band intensity profiles with interactive adjustments for peak detection and background subtraction. Can pick multiple lane areas for faster quantification and analysis.
- **Standard Curve Analysis:** Generate standard curves from known quantities and automatically calculate the amount of protein in sample lanes.
- **Molecular Weight Prediction:** Predict the molecular weight of unknown bands using standard marker sets.
- **Protein Analysis Tools:** Analyze protein sequences for potential N-glycosylation sites and predict the molecular weight of oligomeric or glycosylated forms, with an option to overlay these predictions on your image.
- **Non-destructive Workflow:** Most image adjustments are previews. Permanent changes like cropping, padding, and rotation are undoable.
- **One-Step Auto Analysis:** Draw a single boundary around your gel and let the software automatically straighten the image, detect lanes and bands, and run full densitometry through a guided wizard.
- **Measurement & Calibration Tools:** Calibrate the image to a known length and measure distances and areas directly on the image in real-world units.
- **Flexible Export:** Save your work as a high-resolution image with annotations, a separate configuration file, or export quantitative data directly to Microsoft Excel. Keeps the transparency of the images for posters or other illustrations.

2. Getting Started

2.1 Installation & Dependencies

Gel Blot Analyzer is a Python-based application. There are multiple ways to run the application. You can download the compiled application. However, since the application

is not signed so running it on Windows or MacOS will show an error that this app is from an unknown publisher. As long as you downloaded it from GitHub Releases, the application is tested and safe to use. However, if you want to run the code and have programming skills, you can compile the application or run it directly as well.

Option 1: Download the Pre-built Application (Recommended)

The easiest way to get started is to download the latest pre-built application for your operating system.

1. Go to the [Releases](#) page of this repository.
2. Download the newest release based on numbering:
 - **For macOS:** Download the .zip file.
 - **For Windows:** Download the .zip file.
3. Unzip the file and run the application. No installation is needed.

NOTE: Since the application for MacOS is not signed it may not run or show damaged application, so copy the .app file to the Applications Folder, open terminal and run this command **xattr -cr "/Applications/Gel Blot Analyzer.app"**. For Windows, since it is not from a popular publisher it may show an error that it is blocked by Windows Smart Screen (Windows Defender) so open the application and click on More Info and click on Run Anyway.

Option 2: Run from Source (For Coders)

If you are a coder and want to run the script directly, you can use Conda for a stable environment.

Using Anaconda (Recommended for Stability)

1. Clone the Repository/Download the Library:

```
git clone https://github.com/Anindya-Karmaker/Gel-Blot-Analyzer.git
cd Gel-Blot-Analyzer
```

Or you can download the whole repository as .zip file by clicking on the green button named **<> Code** and select **Download Zip** on the main GitHub Page

2. **Create a New Conda Environment:** This command uses the provided file to create a new environment named **gel_analyzer_env** with all the correct library versions.

```
conda env create -f environment.yml
```

3. Activate the Environment:

```
conda activate gel_analyzer_env
```

4. Run the Program:

```
python Gel_blot_analyzer.py
```

NOTE: If you get any error or crashes particularly due to missing library or C-code header mismatch (usually happens on Windows- Anaconda library issue), run this code:

```
pip install --force-reinstall --no-cache-dir numpy scikit-image
```

Using pip

You can also use pip with the provided requirements.txt file, although dependency conflicts can be more common.

1. Clone the Repository/Download the Library:

```
git clone https://github.com/Anindya-Karmaker/Gel-Blot-Analyzer.git  
cd Gel-Blot-Analyzer
```

Or you can download the whole repository as .zip file by clicking on the green button named **<> Code** and select **Download Zip** on the main GitHub Page.

2. Create and activate a virtual environment (Recommended):

```
python -m venv venv  
# On Windows  
.\venv\Scripts\activate  
# On macOS/Linux  
source venv/bin/activate
```

3. Install requirements:

```
pip install -r requirements.txt
```

4. Run the program:

```
python Gel_blot_analyzer.py
```

NOTE: If you get any error or crashes particularly due to missing library or C-code header mismatch (usually happens on Windows- Anaconda library issue), run this code:

```
pip install --force-reinstall --no-cache-dir numpy scikit-image
```

2.2 Building the Application from Source

If you wish to compile the application yourself, you will need pyinstaller.

1. **Activate your Python environment** (either Conda or venv) where you have installed the dependencies.

2. **Install pyinstaller:**


```
pip install pyinstaller
```

3. **Run the Build Command:** This project includes pre-configured .spec files for both macOS and Windows. Using these is the most reliable way to build.
 - **On macOS:** pyinstaller build_macos.spec
 - **On Windows:** pyinstaller build_windows.spec

The final, self-contained application will be located in the **dist** folder.

2.3 The Main Interface


The main window is divided into three primary sections: Toolbar, Image Viewer, and Image Control/Adjustment Tabs.

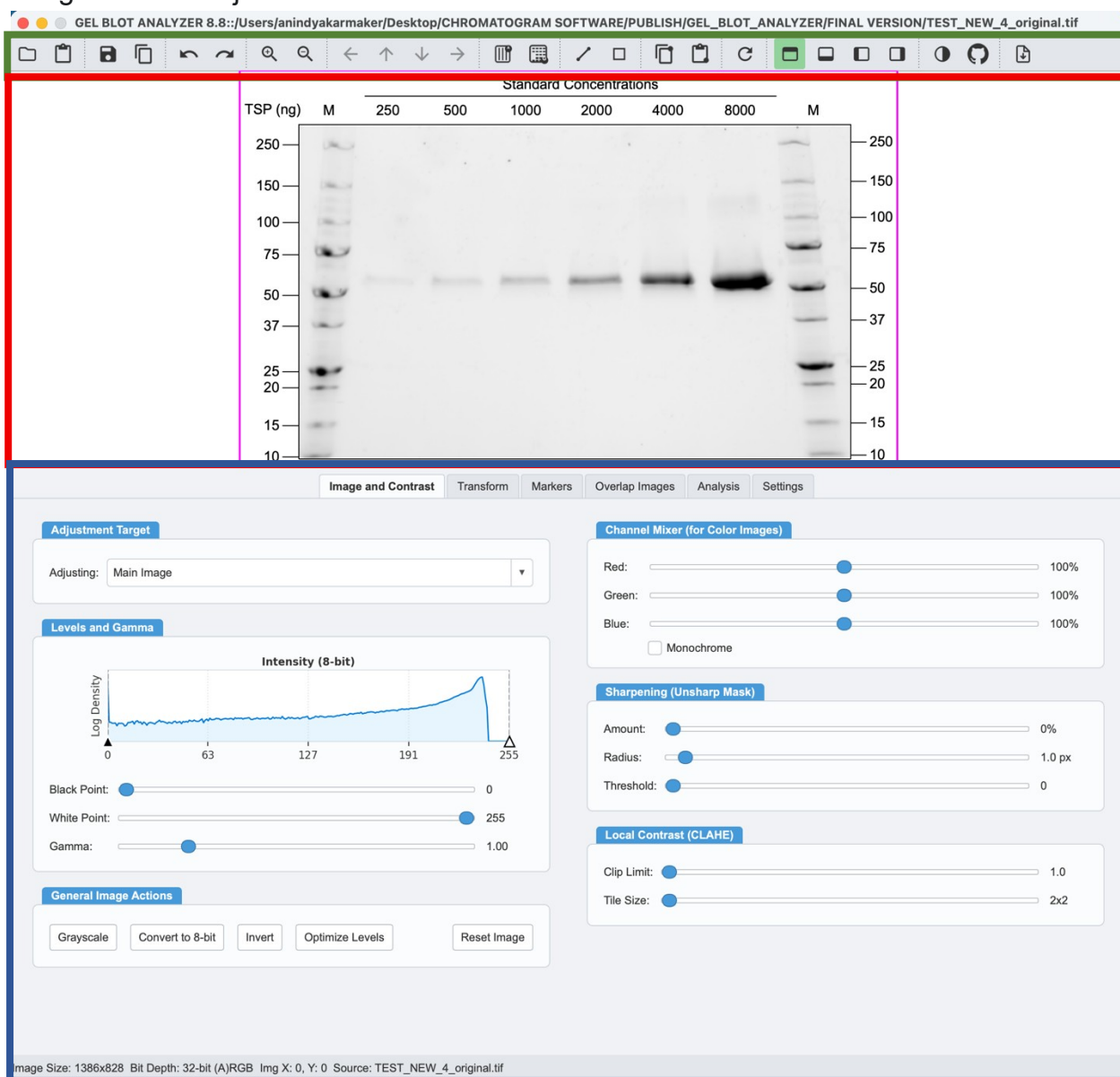


1. **Toolbar:** Provides quick access to the most common functions like loading/saving files, undo/redo, zoom, pan, and annotation tools.

Functions	Icons (v8.0)	Icons (v7.0)
Open, Paste from Clipboard	 	 
Save, Copy to Clipboard	 	 
Undo, Redo	 	 
Zoom In, Zoom Out	 	 
Pan Left, Pan Top, Pan Bottom, Pan Right	   	   
Automatic Lane Marker Labeling		
Auto Gel Analysis		
Create Line, Create Rectangle	 	 
Copy Marker Layout, Paste Marker Layout	 	 
Reset Changes		
Image Viewer Location: Top (Selected), Bottom, Left or Right	   	   
Dark/Light Mode		
Open GitHub Link		
Save as Microsoft® PowerPoint		

2. **Image Viewer:** The central area where your image is displayed. This is an interactive canvas where you will place markers, draw shapes, and define analysis regions. You can change the location of the Image Viewer by clicking on the Image Viewer Location icons. The **magenta** boundary box indicates the rendering region which will be exported when you save image or copy the image to clipboard.

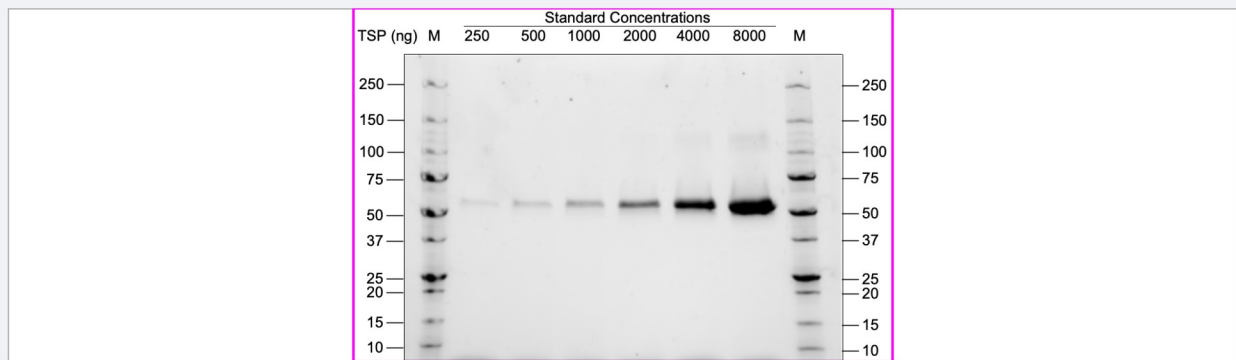

3. **Settings Tabs:** A collection of tabs that organize all the detailed controls for Image and Contrast adjustment, Transform and pad images, Markers placement, Overlap Images (for blending/combining different images), Analysis for densitometric and protein analysis, and a Settings tab for application preferences such as CPU/GPU processing, UI scaling, and window dimensions.

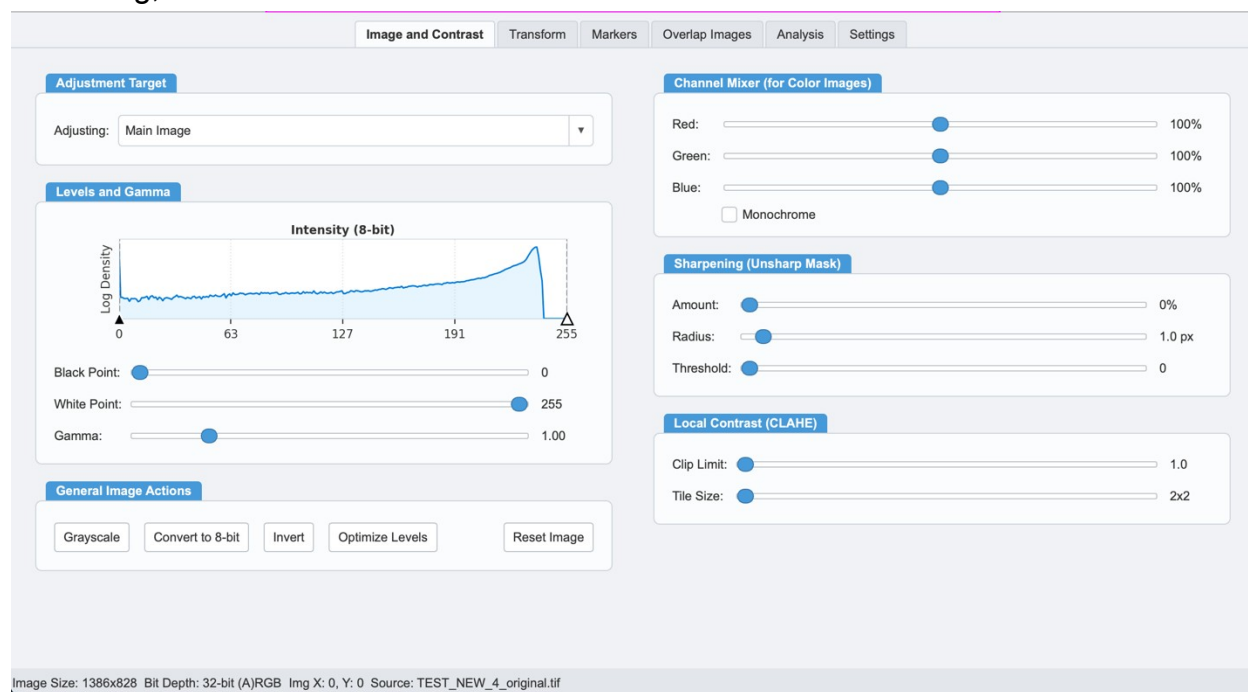

3. Core Concepts

3.1 The Undo/Redo System

Nearly every action that modifies the image or its annotations (e.g., adding a marker, cropping, applying padding, rotating) is saved to an undo stack. It is possible to go back up to 20 steps. Image contrast settings can be reset/adjusted infinitely.

Undo: Use the undo button on the toolbar or the shortcut (Ctrl+Z / Cmd+Z) to reverse the last action.

Redo: Use the redo button or shortcut (Ctrl+Y / Cmd+Shift+Z) to redo the action.

3.1 Saving Your Work

The software uses a paired-file saving system to preserve all your work, including annotations and settings. When you save, it creates three files:

- **[basename]_original.png:** A copy of the modified image *before* any annotations were added.
- **[basename]_modified.png:** A high-resolution image with all your markers and annotations "burned in".
- **[basename]_config.txt:** A text file containing all the data for your markers, shapes, and settings.

When you load the _original.png file, the software will automatically find and apply the corresponding _config.txt file, restoring your session exactly as you left it. Even if you copy the image and paste it, based on the labeling of the file the configuration will be loaded and applied.

Note: Please do not change the file name by removing **_original** or **_config** otherwise the program will fail to load the files properly.

4. Basic Operations: File & View

4.1 Loading an Image

Action: Click the "Load Image" button on the toolbar or use the shortcut (Ctrl+O / Cmd+O).

Function: Opens a dialog to select an image file (.png, .jpg, .tif, etc.). If a corresponding _config.txt file is found, it will be loaded automatically. You can load even raw 16-bit images but the raw images may need further adjustments like image invert operation and adjusting the black and white point values in the **Image and Contrast** tab.

4.2 Saving an Image

Action: Click the "Save with Config" button or use the shortcut (Ctrl+S / Cmd+S).

Function: Opens a dialog to save your work using the paired-file system described in [Core Concepts](#). It is possible to save the image using different formats (PNG, JPG, TIFF, etc). PNG is ideal for pasting into presentations or documents where transparency is necessary and TIFF is ideal for lossless images ideal for publications.

Save to PowerPoint (.pptx): Use **Save to PPTX...** (on the toolbar and in the **File** menu) to export the figure as an editable PowerPoint file. The gel/blot image goes in as a single high-resolution background picture, and every annotation — the Left/Right/Top markers, custom text markers, and custom shapes (lines, rectangles and arrows) — is added as a native, fully editable PowerPoint text box or shape. You can click any label in PowerPoint to change its text, font, size or colour (or copy the whole figure from PowerPoint into Word). This is ideal when you need to tweak annotations after export without re-rendering in GBA.

4.3 Copy & Paste

Copy Image (Ctrl+C / Cmd+C): Copies the currently displayed image *with all visible annotations* to the clipboard as a high-resolution PNG in the clipboard for pasting directly into other software like Microsoft® Word or PowerPoint etc.

Paste Image (Ctrl+V / Cmd+V): Pastes an image from the clipboard into the viewer, replacing the current image. If the clipboard contains a path to an image file, it will load that file and its associated config if available.

4.4 Navigating the Image: Zoom & Pan

Zoom In/Out:

Use the + and - buttons on the toolbar.

Use the mouse wheel or 2-finger pinch on a touchpad to zoom in or out the image.


Use shortcuts (**Ctrl** + = or **Cmd** + = for zoom in, **Ctrl** + - / **Cmd** + - for zoom out).

Pan (Move the View):

When zoomed in, **click and drag** the image with the **right mouse** button.
Use the arrow buttons on the toolbar or use the arrow keys on your keyboard.

4.5 Changing the Layout

You can change the position of the Image Viewer relative to the Settings Tabs for your convenience. Use the layout buttons on the toolbar to position the viewer on the Top, Bottom, Left, or Right.

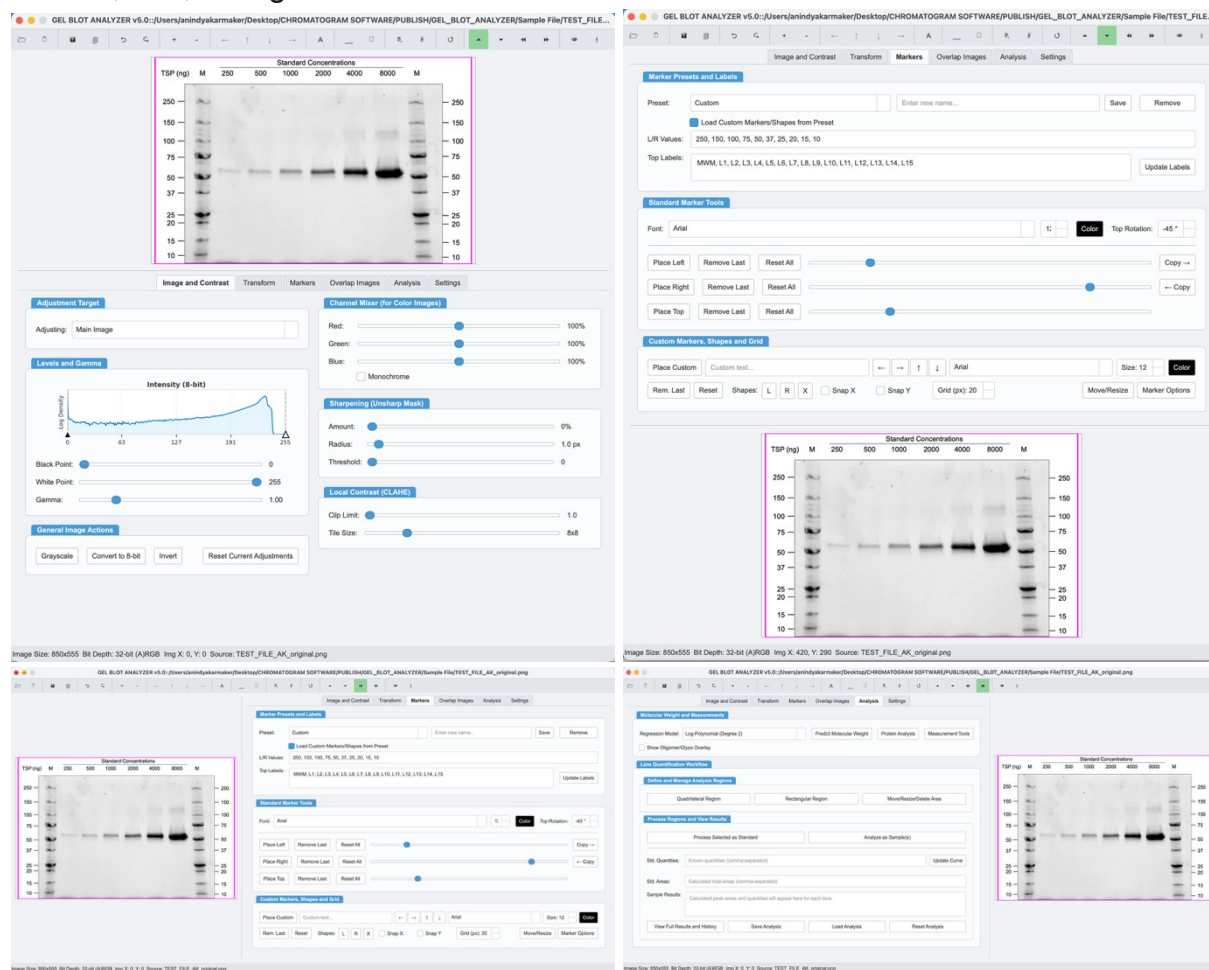

4.6 Adjusting the User Interface

You can control the size of the Image Viewer Area (Viewer Fixed Width/Viewer Fixed Height) and the width of the control section using the Program UI width. Moreover, if the icons/symbols appear too or small, please adjust the UI scaling from 0% to 250% from the Settings Tab. Note: Some settings require restart of the application to work properly. The application by default uses GPU for acceleration but it can be toggled off to use mostly CPU.

Image and Contrast

Transform

Markers

Overlap Images

Analysis

Settings

Processing Unit Priority

☒ GPU (Graphics Processing Unit)

☐ CPU (Central Processing Unit)

Target GPU ID:

:GPU:0

(Currently Active: GPU: Apple M1 Pro)

User Interface Scaling (Requires Restart)

UI Scale:

100% (No Scaling)

Application Dimensions (Requires Restart)

Viewer Fixed Width:

550

Viewer Fixed Height:

350

Program UI Width:

1000

Apply and Save Settings

Reset All Settings

image Size: 1514x813 Bit Depth: 32-bit (A)RGB Img X: 0, Y: 0 Source: TEST_NEW_3_original.tif

14

5. Image Adjustments & Transformations

5.1 Image and Contrast tab


The **Image and Contrast** tab provides a suite of tools to optimize the visual appearance of the image for analysis. With the exception of "Grayscale", "Convert to 8-bit" and "Invert," all adjustments are non-destructive previews, meaning they do not permanently alter the original image data until you perform an action like rasterizing.

Adjustment Target: This dropdown menu is the most important control on this tab. It selects which image is currently being adjusted: the **Main Image**, or one of the overlays (**Overlay 1** or **Overlay 2**). This allows you to apply different contrast, color, and sharpening settings to each image independently, which is essential for matching the appearance of different images (e.g., a gel and a western blot) when creating a composite overlay (for more details see **Overlap Images** tab).

Levels and Gamma: These sliders control the overall brightness and contrast.

Black/White Point: Adjusts the dynamic range. Pixels darker than the **Black Point** become pure black, and pixels lighter than the **White Point** become pure white. This is the primary tool for maximizing contrast. The effect is visualized by the markers on the **Log-Histogram** window.

Gamma: Adjusts the brightness of the image's mid-tones without affecting the absolute black and white points, helping to reveal faint details.

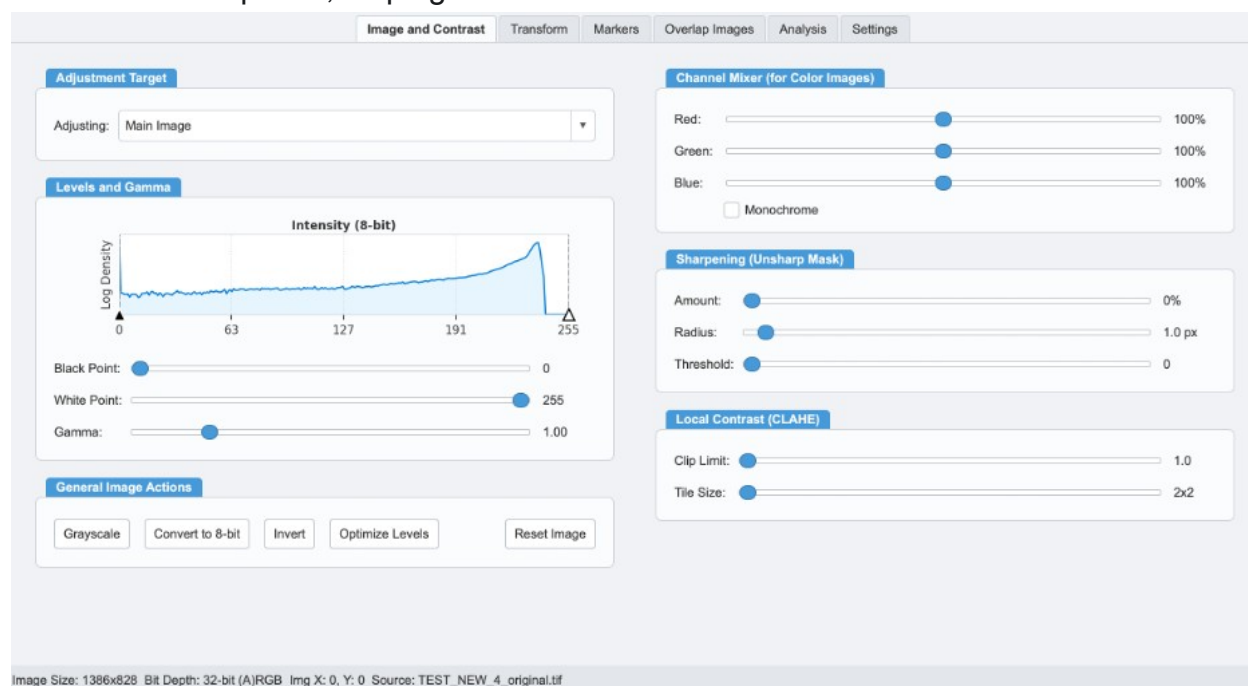

Channel Mixer (for Color Images): This tool gives you precise control over the color channels in an RGB image. Increase or decrease the contribution of each color channel (RGB). This is useful for isolating signals from specific fluorescent channels or for creating a custom color balance. **Monochrome:** Creates a grayscale image based on the weighted contributions of the red, green, and blue channels that have been set.

Sharpening (Unsharp Mask): Increases the contrast along the edges of features, making them appear sharper. **Amount:** Controls the strength of the sharpening effect. **Radius:** Determines the size of the area around the edges that will be affected. **Threshold:** Prevents sharpening in low-contrast areas, which helps to avoid amplifying image noise.

Local Contrast (CLAHE): Contrast Limited Adaptive Histogram Equalization improves contrast by enhancing it locally across different regions of the image, rather than globally. This is excellent for observing or visualizing very faint bands without causing saturation of other bands. **Clip Limit:** Controls how much contrast is amplified, preventing noise from being excessively brightened. **Tile Size:** Defines the size of the local regions for histogram equalization. Higher is better.

General Image Actions:

Grayscale (Ctrl+B / Cmd+B): A permanent action that converts the image to a high-precision grayscale format, which is recommended for most quantitative analyses. Can use Undo option to revert back the changes.

Convert to 8-bit: If your program is running really slow and you do not need increased precision, this option can lower the image resolution for easier operation. However, quality may be affected. Can use Undo option to revert back the changes.

Invert (Ctrl+I / Cmd+I): A permanent action that inverts all image colors (e.g., black to white). This is essential for densitometry specially for DNA/RNA gel analysis, as it transforms dark bands into positive peaks for analysis. Can use Undo option to revert back the changes.


Optimize Levels: Auto adjusts the black and white point for optimum display and analysis. May not work with some images with multiple peaks.

Reset Adjustments: Resets all preview sliders on this tab to their default values for the currently selected adjustment target.

NOTE: Densitometric analysis applying most image or contrast options particularly **Channel Mixer** or **CLAHE** may change the band analysis significantly and user caution is advised.

5.2 Transform tab

The **Transform** tab contains all the tools needed to geometrically correct and prepare your image for analysis. It is organized into two main sections: non-destructive **Global Transformations** that can be viewed as live previews and then applied, and destructive **Cropping and Padding** operations that permanently alter the image canvas (can undo changes).

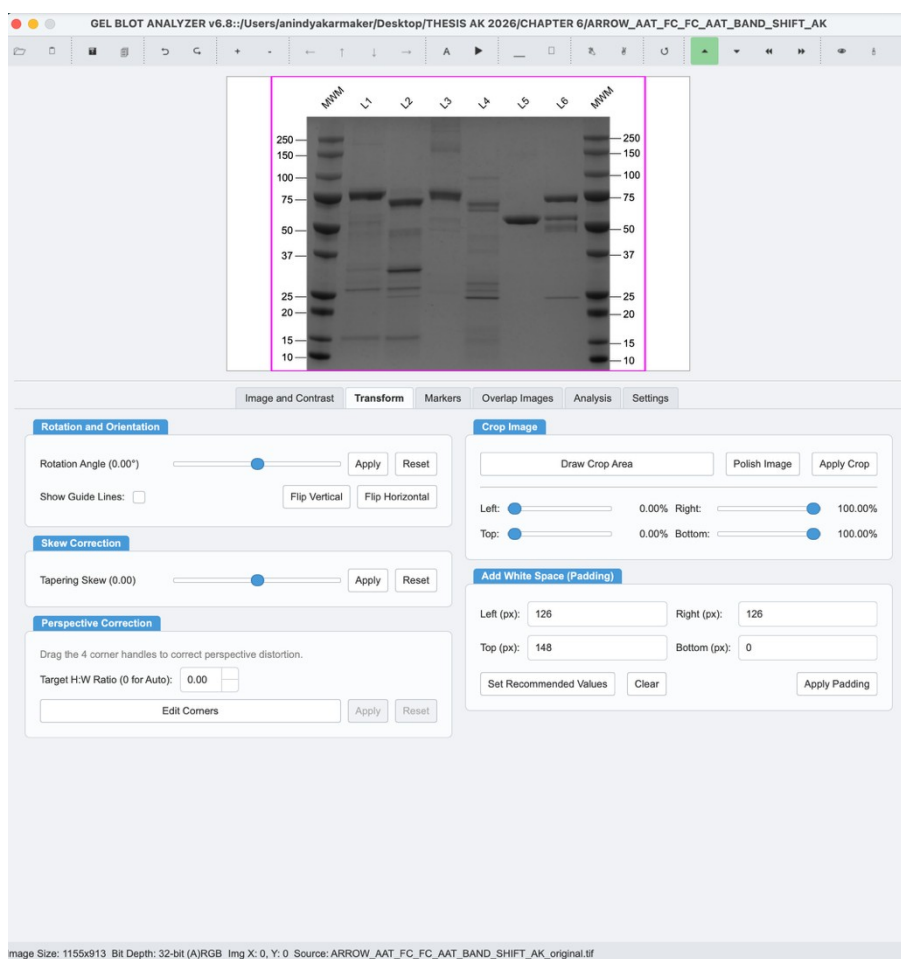

5.3 Global Transformations (Live Previews)

These tools allow you to see the effect of your changes in real-time before applying them.

Rotation and Orientation: Use the **Rotation Angle** slider for a live preview to straighten a tilted gel or blot. You can toggle the **Show Guide Lines** checkbox (Shortcut: Ctrl+G) to display a crosshair, making precise alignment easier. Once you are satisfied, click **Apply** to make the rotation permanent. Note that applying a rotation will remove any standard L/R/Top markers, as their alignment becomes ambiguous. The **Flip** buttons provide immediate, permanent vertical or horizontal mirroring of the image.

Skew Correction: This tool corrects for perspective distortion, often seen as "smiling" or "frowning" bands in a gel. Use the **Skew** slider for a live preview of the correction. When the lanes appear parallel, click **Apply** to make the transformation permanent.

5.4 Destructive Edits (Permanent Actions)

These actions modify the fundamental dimensions of your image. It is recommended to perform these edits before adding detailed annotations.

Crop Image: This is used to remove unwanted areas of the image.

- **Draw Crop Area:** The primary method. Click this button, then click and drag on the image to define the area you want to keep.
- **Sliders:** For precise adjustments, you can use the sliders to define the crop boundaries as a percentage of the image dimensions.
- **Apply Crop:** This is the final step. Clicking this button permanently removes everything outside the defined crop area. All existing markers and shapes will be repositioned relative to the new top-left corner.

Add White Space (Padding): This permanently adds a blank border around your image. This is extremely useful for creating clean space for markers and other annotations, ensuring they don't cover important image data. You will notice the changes as the magenta boundary region will change highlighting the change in the rendered output.

- **Input Fields:** Enter the number of pixels to add to each side.
- **Set Recommended Values:** A convenient button that suggests default padding values based on the current image size.
- **Apply Padding:** Click this to finalize the action. The image canvas will be enlarged, and all existing markers and shapes will be automatically shifted to maintain their correct positions relative to the original image content.

Perspective Correction: The user picks four corner points based on the level of skewness of the gel/image and it can auto-adjust the final image into a perfect rectangular shape with the perspective corrected.


- **Target H:W ratio:** Default is set to auto (0) but the user can make it can fix height to width ratio, useful if you know the gel height to width ratio.

Note: Perspective correction and/or skew correction modifies the bands significantly. This feature is significantly useful to get accurate molecular weight marker prediction, but for densitometry, it can change the values. User caution advised.

6. Annotations: Markers & Shapes

6.1 Markers tab

The **Markers** tab is the control center for all annotations.

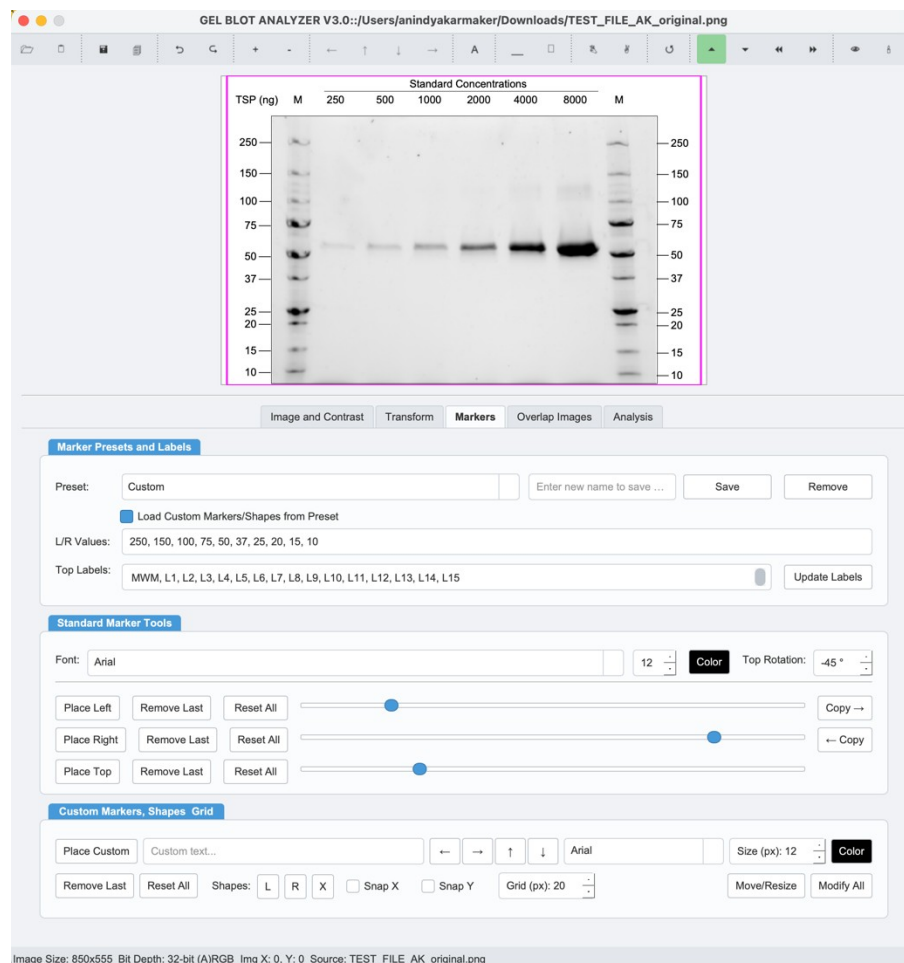

Standard Markers (Left, Right, Top)

These are typically used for molecular weight standards or lane labels.

1. **Select a Preset:** Choose a standard from the **Preset** dropdown (e.g., a Bio-Rad or Thermo ladder). This populates the internal list of values. If you want to create your own markers then choose the **Preset** that is closest to your target and then choose the **Preset: Custom** and enter or edit the values separated by comma (,) and give a suitable name and click on **Save Preset**. The save preset function also saves the **Top Marker Labels** as well. The **Save Preset** will also save any **Custom Markers** as well. If you load a **Preset** but do not want to load the saved **Custom Markers**, please uncheck the option **Load Custom Markers/Shapes from Preset**.
2. **Place Markers:**
 - Click **Place Left** (Shortcut: Ctrl+Shift+L / Cmd+Shift+L).

- Click on the image at the vertical position of each corresponding marker band, from top to bottom. The software will automatically apply the labels from the preset.
 - Repeat for **Place Right** (Ctrl+Shift+R / Cmd+Shift+R) and **Place Top** (Ctrl+Shift+T / Cmd+Shift+T).
3. **Copy Markers:** If you want to duplicate the markers from left side or right side, click the (Copy →) or (Copy ←) button.
 4. **Adjust Offset:** Use the **Offset** sliders to move all markers on one side horizontally (for L/R) or vertically (for Top). Sometimes duration image manipulation if the markers end up going outside the preview window then use it to correct the marker locations.
 5. **Remove/Reset:** Use the **Remove Last** to remove last marker entry or **Reset All** buttons to reset the markers.

You can also modify the font type, size, color and rotation (for top labels) by changing the variables in this tab.

NOTE: If **Line X** or **Line Y** is on along with the **Snap**, the markers and other parameters will be constrained along the center of the grid lines.

6.2 Custom Markers & Shapes

Use this section for any other annotations like placing arrows or molecular weight, etc.

1. Custom Text/Arrows:

- Enter text in the text box or click an arrow button from the default presets.
- Select the desired font, size, and color.
- Click Place Custom.
- Click on the image to place the marker.

Note: You can change the color, text, size and font of the custom markers and can see an overlay before placing it. Please click **Place Custom** if you change the parameters to see the updated overlay of the custom text/shape.

2. Shapes (Lines, Rectangles & Arrows):

- Select the desired size (thickness) and color.
- Click the **L** (Line), **R** (Rectangle) or **A** (Arrow) button in the “Shapes” section.
- Click and drag on the image to draw the shape. Hold **Shift** while drawing to constrain lines to 45° angles and rectangles to squares.
- **Arrows:** The **A** tool draws an arrow in any direction — the arrowhead points to wherever you release the mouse. Hold **Shift** to constrain the angle to 0°, 45° or 90°. In **Move/Resize** mode, drag the two end handles to re-aim or lengthen the arrow, and drag the two side handles (perpendicular to the shaft) to make it thicker or thinner.

3. **Move/Resize Custom Items:** Click the **Move/Resize** button to enter interaction mode. You can then click and drag markers or shapes to reposition them. For shapes, drag the corner handles to resize.

NOTE: If **Line X** or **Line Y** is on along with the **Snap**, the markers/shapes will be constrained along the center of the grid lines.

6.3 Managing Custom Markers/Annotations

You can use the **Marker Options** to open the **Modify Custom Markers and Shapes** window to adjust parameters like Text/Label, Location, Font Type, Font Size and other properties.

	Type	Text/Label	Coordinates	Font	Size	Bold	Italic	Color	Actions
1	Marker	M	127.4,42.6	Arial	12	☐	☐	#000000	Delete
2	Marker	M	703.5,42.6	Arial	12	☐	☐	#000000	Delete
3	Marker	TSP (ng)	52.8,42.6	Arial	12	☐	☐	#000000	Delete
4	Marker	250	194.1,42.6	Arial	12	☐	☐	#000000	Delete
5	Marker	500	283.0,42.6	Arial	12	☐	☐	#000000	Delete
6	Marker	1000	363.9,42.6	Arial	12	☐	☐	#000000	Delete
7	Marker	2000	451.2,42.6	Arial	12	☐	☐	#000000	Delete
8	Marker	4000	535.3,42.6	Arial	12	☐	☐	#000000	Delete
9	Marker	8000	624.2,42.6	Arial	12	☐	☐	#000000	Delete
10	Marker	Standard Concentrations	433.7,14.1	Arial	12	☐	☐	#000000	Delete

The **Shift X** and **Shift Y** shift the overall markers X and Y axis respectively. The **Scale X** and **Scale Y** shrink or expand the custom markers. The **Font Scale** adjusts all of the custom marker fonts globally except the Top/Left/Right markers since they have their own font/color and other properties settings.

Copy/Paste Custom Items (Experimental Feature): Use the copy/paste buttons on the toolbar to copy all custom markers and shapes from one image, load the other image and then paste them onto another. Use the Shift X/Y to change the location of the pasted markers and Scale X/Y to control their size. Font Scale changes the global font size for everything.

Note: Only use the copy/paste after placing the left/right/top markers otherwise the custom markers will be reset.


6.4 The Grid System & Snapping

- **Show Grid:** Check **Line X** (Shortcut: Ctrl+Shift+X or Cmd+Shift+X) and **Line Y** (Shortcut: Ctrl+Shift+Y or Cmd+Shift+Y) to display a grid.

- **Grid Size:** Adjust the grid spacing in pixels using the **Grid (px)** spinbox (Shortcut: Ctrl+Shift+G or Cmd+Shift+G to enable the grid). To increase or decrease the grid size, use the **Grid (px)** spinbox (Shortcut: Ctrl+Shift+ ↑ or Cmd+Shift+ ↑ and Ctrl+Shift+ ↓ or Cmd+Shift+ ↓)
- **Snap:** When enabled (Shortcut: Ctrl+Shift+N or Cmd+Shift+N), every mouse click — placing markers, drawing shapes, etc. — snaps to the nearest grid intersection. Turning Snap on automatically shows the grid so you can see what you are snapping to: if neither Line X nor Line Y is currently displayed, both are switched on; if one is already on, your choice is kept and snapping simply uses it. Turning Snap off leaves the grid lines as they are.

7. Overlap Images

You can overlap a WB on top of the SDS-PAGE or combine multiple images using this function.

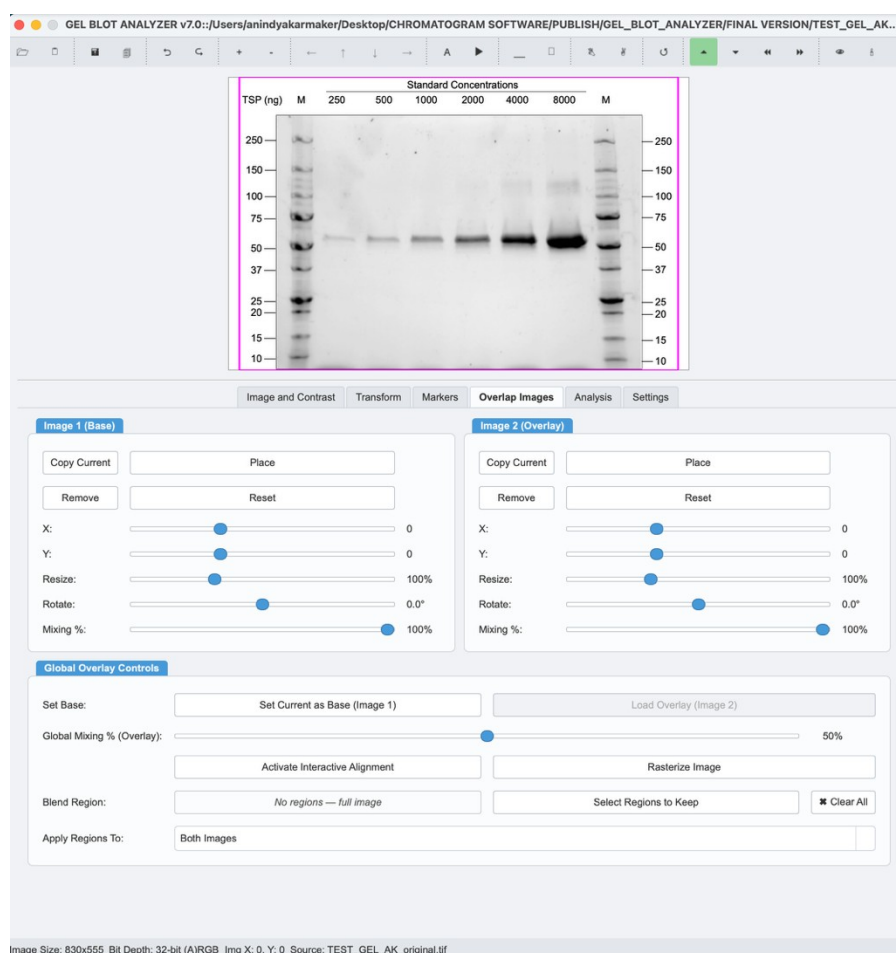

1. Paste an image into the image viewer and make any edits using the previous tabs. Once you are happy with the changes, click on **Copy Current Image** in the **Image 1 (Base)** section and it will copy it into memory.

2. Load the second image if you have a second image, make the changes and click on **Copy Current Image** in the **Image 2 (Overlay)** section.
3. Now load the final or last image on top of which you want to place the other two images. Use the **Padding** tab to create the required spacing to place the other two images or if you want to overlap then on top of each other.
4. Click on **Place** in the **Image 1 (Base)** section and/or **Place** in the **Image 2 (Overlay)** section. Place adds the layer on top of the image you currently have open — it does not replace it — so you can build up a figure from several gels and blots. (The one exception is **Set Current Image as Base** below: because that makes Image 1 a copy of the image already on screen, the copy is shown instead of the original rather than on top of it, so it is not drawn twice.)
5. You can use the sliders to adjust the size, position, rotation or level of blending. You can click on **Activate Interactive Overlay Mode** and you can use it to resize or move the image by dragging with mouse.
6. **Remove** takes a layer off the canvas but keeps it loaded, so you can reposition or re-stack it and press **Place** again without copying or loading it a second time. Use **Reset** when you actually want to empty the Image 1 or Image 2 slot — **Reset** is the only control that discards the stored image.
7. Once you are happy with the changes, please click on **Rasterize Image** to finalize it. Then you can go ahead and add the left/right/top markers or add more padding as required. Rasterizing takes the placed layers off the canvas (they are now part of the image itself) but leaves them loaded in their slots, so you can immediately **Place** the next panel against what you just baked in and rasterize again — repeat this to assemble a multi-panel figure. **NOTE:** Rasterizing bakes the layers permanently into the image. It is a single undo step, so Ctrl+Z (Cmd+Z) puts the layers back exactly as they were placed — but once you save or move on, the merge cannot be taken apart again.
8. If you just want to merge two images like Chemiluminescence and Colorimetric images for WB (does not need individual adjustments), you can load the first image and click on **Set Current Image as Base (Image 1)** and then click on **Load Overlay (Image 2)** to load the second image.
9. You can use **Blend Region** to create selection rectangle to keep areas that you want to blend/overlap and use the **Apply Regions To** select which areas need to be blended.

Note: Default blending is set to 50% for global and 100% for individual. If you want to overlap images on-top this behavior is perfect. However, if you want overlap one image in an empty spot on another image, please set the blending to 100%.

Overlap Images — what each control does

Each of the two panels, **Image 1 (Base)** and **Image 2 (Overlay)**, has the same set of controls:

Copy Current — copies the image currently open into that slot's buffer. It does not put anything on the canvas yet; the image is simply held, ready to Place. The snapshot is taken as the image is, so a Rotation Angle you have dialled in but not yet applied is NOT baked into it.

Place — puts the buffered image onto the canvas at the slot's X/Y position. It is added on top of the image you currently have open rather than replacing it, so you can build a figure from several gels and blots. (Set Current as Base is the exception — see below.) A layer is clipped to the image's own edges on all four sides: anything you push beyond them is not shown, because Rasterize and Save would discard it. Add padding in the Transform tab to make room first.

Remove — takes the layer off the canvas but keeps it loaded in the slot, so you can reposition or re-stack it and press Place again without copying or loading it a second time.

Reset — empties the slot completely, discarding the stored image. This is the only control that throws the image away; use Remove if you only want it off the canvas.

X / Y — the layer's position in image pixels. Negative values move it off the left/top edge.

Resize — scales the layer. Arrow keys nudge it in 0.01% steps once the slider has focus, which is how you match a layer to the image underneath.

Rotate — turns that layer alone. Rotating the main image in the Transform tab does not affect placed layers; this is the control for turning one.

Mixing % — that layer's own opacity. Image 1 uses this value alone; Image 2 is multiplied by the Global Mixing % as well.

The **Global Overlay Controls** below them apply to the composition as a whole:

Set Current as Base (Image 1) — copies the current image into Image 1 and places it in one step. Because Image 1 is then a duplicate of the image already open, it is shown INSTEAD of it rather than on top of it, so the same picture is not composited twice.

Load Overlay (Image 2) — loads a second image from disk into Image 2. The file dialog opens in the folder the current image came from.

Global Mixing % (Overlay) — applies to Image 2 only, on top of its own Mixing %. Image 1 ignores it. The default is 50%, which suits overlaying two exposures of the same blot; set it to 100% when you are placing a separate panel into empty space.

Activate Interactive Alignment — click a layer to select it, drag to move, mouse wheel to resize, arrow keys to nudge. Press Esc to leave the mode.

Rasterize Image — bakes the placed layers permanently into the image. The layers come off the canvas (they are part of the image now) but stay loaded in their slots, so you can immediately Place the next panel against what you just baked in and rasterize again. This is a single undo step.


Blend Region / Select Regions to Keep — draw rectangles to restrict blending to chosen areas; Clear All removes them.

Apply Regions To — chooses which layer those regions apply to.

8. Quantitative Analysis Workflow


8.1 Defining Analysis Regions

Go to the **Analysis** tab. You must first define the area(s) you want to analyze.

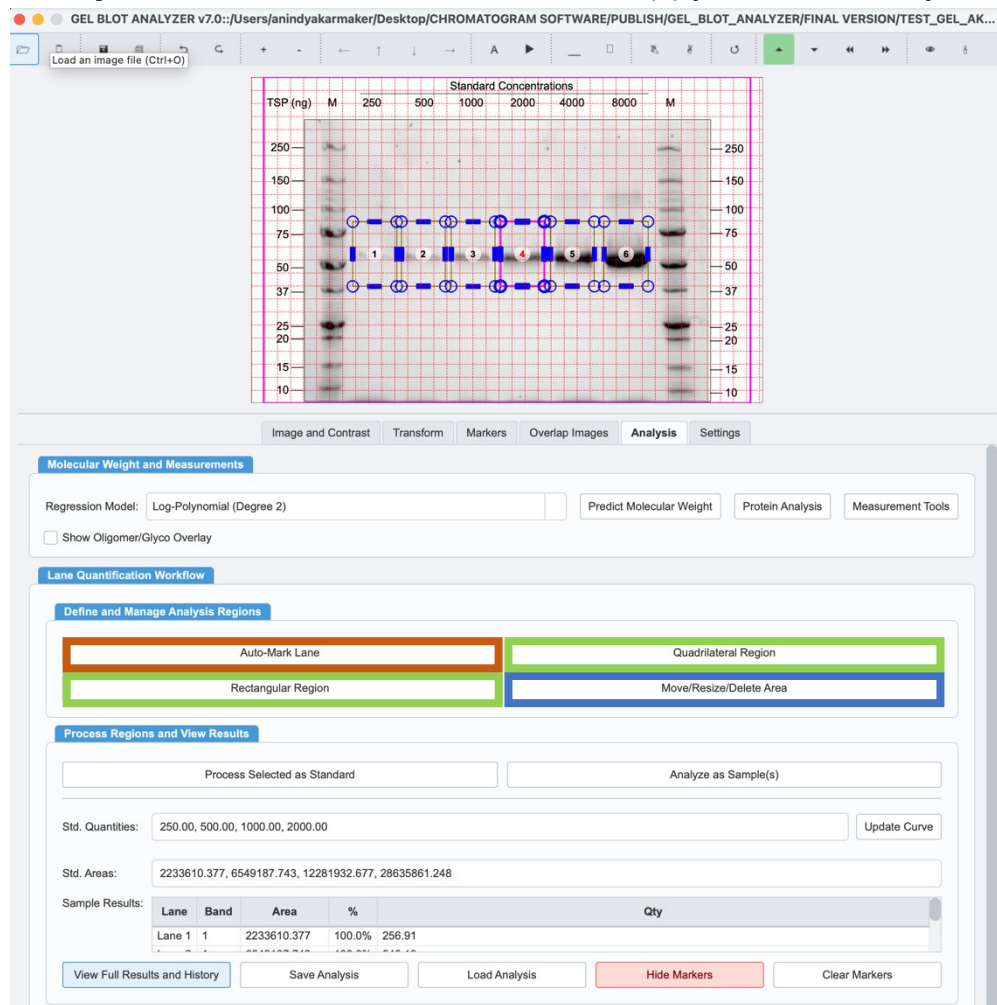

- For picking region/Lane for analysis:
 - Click **Quadrilateral Region** for skewed lanes. Click four points around your lane starting with top left, top right, bottom right and then final bottom left. **Note:** It is important to follow this direction for proper custom area selection.
 - Click **Rectangular Region** for straight lanes. Click and drag to draw a box around your lane.
 - If you want to use the same size of the bounding box, click on **Auto-Mark Lane** and click on the center of the region around which you want the box. It is a much faster way to create the bounding boxes
- Moving/Resizing Regions:** Click **Move/Resize/Delete Area** to enter a mode where you can click a defined lane to select it, then drag it to move or drag its corners to resize.


NOTE: You can use the Snap Grid Functionality to pick uniform areas for analysis since your selection of area impacts the analysis (**CTRL/CMD + Shift + S**). You can use **Shift key + Left Click** Mouse drag to move in straight lines by clicking on the center and you can move edges by clicking on center blue knobs with the same shortcut. You can duplicate the regions by clicking **Ctrl/CMD + Shift + Left Click** Mouse drag from the

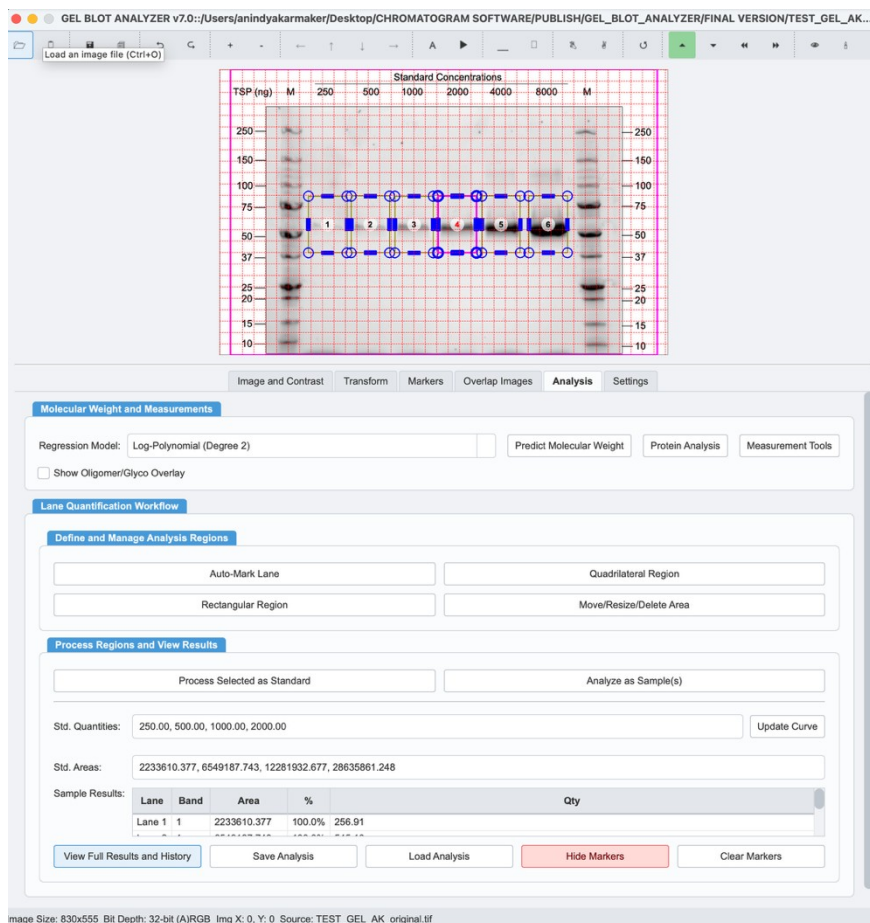

center.

8.2 Processing a Standard:

1. Define a region for your band as discussed before.
2. Use the **Move/Resize/Delete Area** button to highlight the region.
3. Click **Process Selected as Standard**.

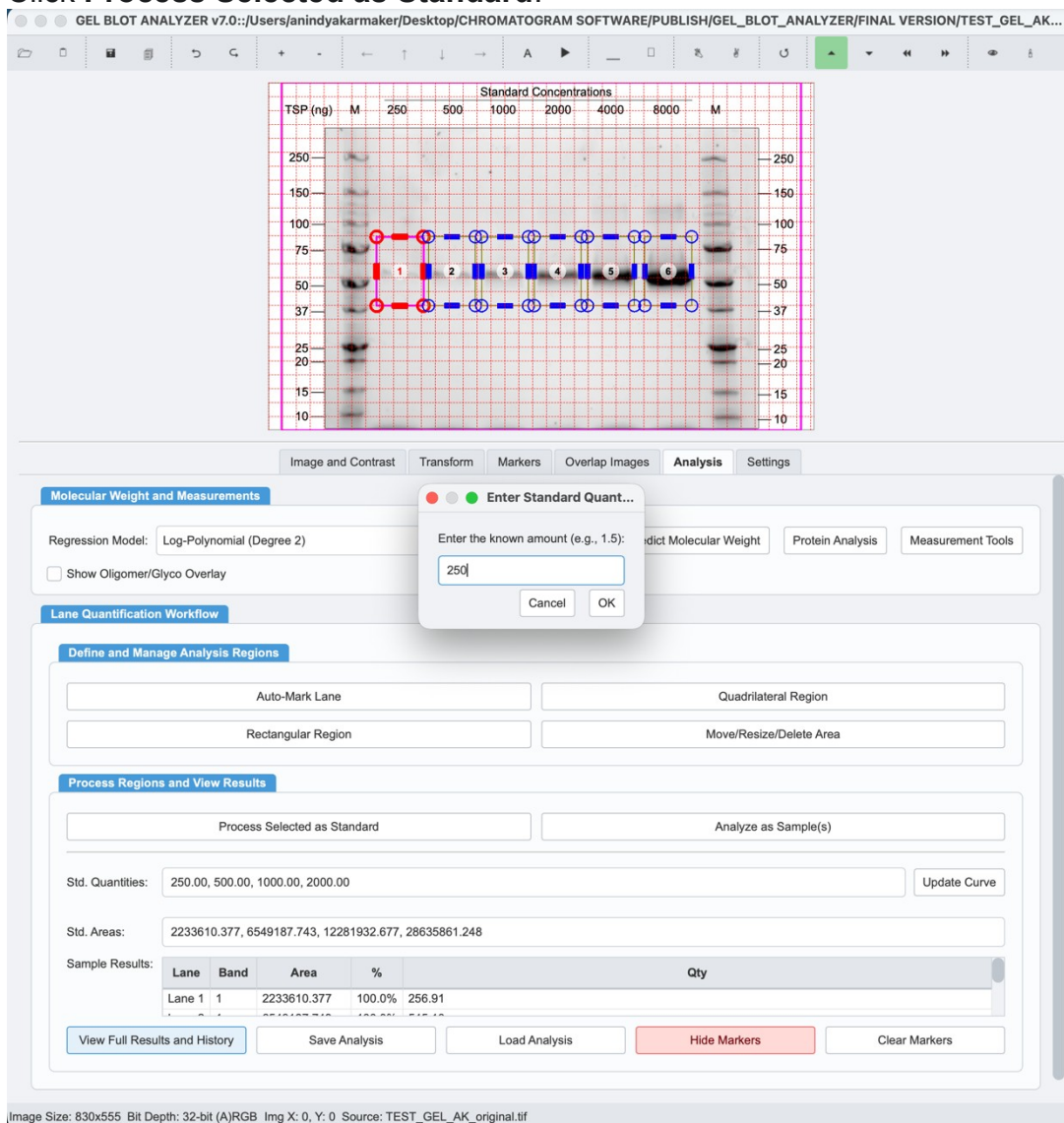


4. A dialog box opens. Enter the known amount of standard and click **OK**. The Adjust Peak Regions and Calculate Areas dialog box opens. It should detect a single peak. Adjust parameters and click **OK**
 5. Repeat for at least two standard lanes to generate a standard curve. The values will appear in the "Std. Quantities" and "Std. Areas" fields. If you want to eliminate custom values, please remove the areas and quantities and click the button **Update Curve**.
- Note:** To get the highest precision in quantification, please make sure to use the same parameters for detection particularly the **rolling ball radius**. Also, the quantification is accurately tested up to 4 μ g of protein loaded onto the gel.

8.3 Analyzing a Sample:

1. Define the region(s) for your sample lane(s).
2. Click Analyze as Sample(s).
3. The Peak Area Dialog will open for each defined lane. Adjust and click **OK**.
4. The software uses the standard curve to calculate the quantity for each band. The results (areas and calculated quantities) will appear in the "Sample Areas" text box.

8.4 Peak Area Dialog

This is where you fine-tune the analysis.

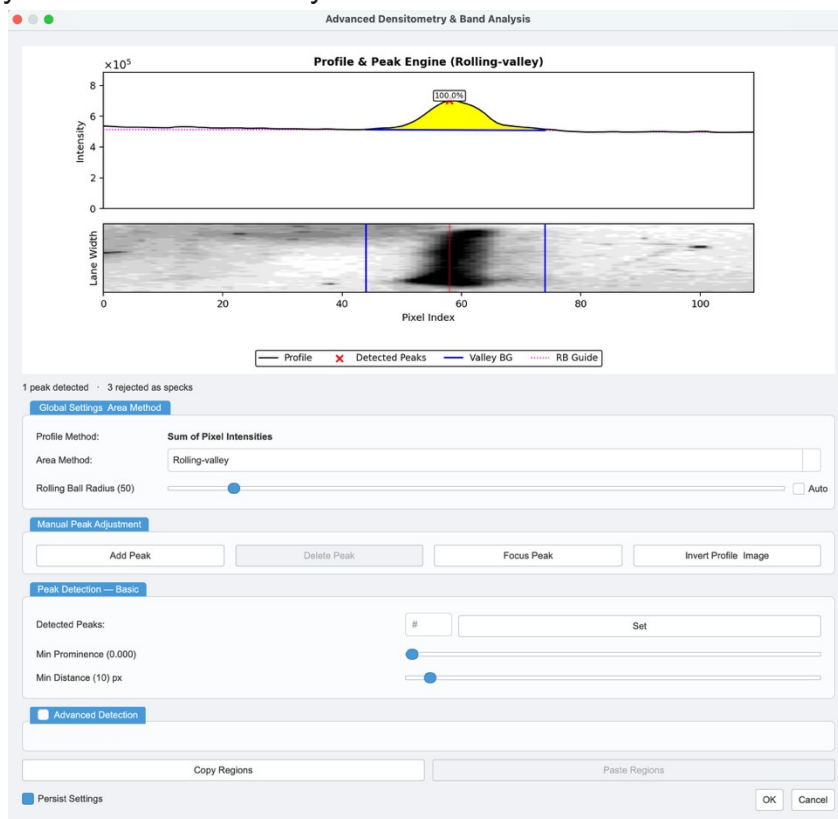

The Profile Plot: Shows the intensity profile along the length of the lane. Detected peaks are marked with 'x', and the calculated area is shaded.

The Image Strip: A preview of your lane. You can directly click and drag the blue boundary lines to manually adjust the integration area for each peak.

8.5 Global Settings Area Method:

Area Method: Choose the background subtraction method:

Rolling Ball: Best for uneven backgrounds. Adjust the **Radius** slider for background subtraction. Best to select **Auto** for best accuracy.

Valley-to-Valley: Defines background as the line between the lowest points (valleys) on either side of a peak.

Rolling-Valley: Defines background as the straight line connecting the start and end of a peak's integration boundaries using rolling ball for background subtraction.

Gaussian Deconvolution: Experimental feature. Good for great resolution SDS-PAGE images.

8.6 Manual Peak Adjustment

Use the **Add Peak** and **Delete Peak** buttons to manually correct the detected peaks directly on the plot. **Focus Peak** is used to highlight the peak relevant boundaries for easier manipulation. **Invert Profile Image** is to invert the peaks to troughs and vice-versa for inverted images

8.7 Peak Detection Settings:

Use the sliders (**Min Prominence**, **Min Distance**) to adjust the automatic peak detection algorithm.

Use the **Advanced Detection** feature to tweak the algorithm for peak detection and denoising. By default, the software does not perform denoising and smoothing which may be required for some analysis and that option is there. Good typical values are **0.5** for denoising and **2.0** for smoothing sigma. Use the Min Height and Min Band Width to remove noise or speck from detection. Good starting values are **0.10** and **6 px** respectively.

Detection control reference (tune roughly in this order): the controls below are split between the Basic group and the collapsible Advanced Detection group.

Min Prominence — how far a peak must rise above the dips on either side, as a fraction of the full trace height. This is the single most effective noise control: raise it first if specks are being counted as bands; lower it to recover genuine faint bands.

Min Distance (px) — the minimum spacing allowed between two peaks. Raise it to stop one broad band being split in two; lower it to resolve closely spaced bands.

Denoise Sigma — blurs the 2-D image (and widens an impulse-removing median filter) before the trace is built. Raise it for grainy/speckled images; too high merges neighbouring bands. A small median filter always runs, even at 0. Typical: 0.5.

Smoothing Sigma — smooths the 1-D trace only. Helps detection on noisy lanes; too high flattens shoulders and can shift peak centres. Typical: 2.0.

Min Height — the minimum peak height, as a fraction of the trace range. Rejects faint bands outright; Min Prominence is usually the better speck filter. Typical: 0.10.

Min Band Width (px) — the minimum peak width at half-maximum. Real bands span several pixels whereas single-pixel spikes do not, so raise this to drop sharp noise spikes in the 1-D trace. Typical: 6 px.

Speck Removal (Advanced Detection): this is the most direct control for the common problem of dust, fibres or staining specks being picked up as bands. Unlike **Min Height** and **Min Band Width** (which act on the 1-D trace), **Speck Removal** works on the 2-D lane image before the trace is built: any bright blob smaller than the chosen area (in px²)

is erased, so it can never become a peak, while genuine bands — which always span many pixels — are left untouched. Leave it at **0** to disable; start around **10–40 px²** and raise it until the specks disappear.

You can also manually set the number of peaks if the detection algorithm fails by changing the number of detected peaks and click on Set and then narrow the region of peak peak dragging the blue boundary lines on the **Image Strip**. But the easier option is to just use **Add Peak** to click on the peaks manually. You can also delete the peaks by clicking on it and press delete key or click on the **Delete Peak** button. The **Focus Peak** is useful if you have too many peaks and cannot bind the boundary edges for area estimation. The **Invert Image Profile** is necessary if your bands are white on a dark background (typical for DNA/RNA gels)

Tip — a good starting point: Rolling-valley background with Auto radius; Min Prominence ≈ 0.05 ; Min Distance ≈ 10 px; Denoise 0.5; Smoothing 2.0; Min Height 0.10; Min Band Width 6 px; Speck Removal 0 (raise to 10–40 px² only if specks are still detected). Tune one control at a time and watch the Profile Plot update live; settings you set on one lane are shared with the others, and you can fix individual lanes afterwards.

NOTE: The band detection is tweaked to work only with dark bands on a white background. Please **Invert** the image in the Image and Contrast Setting or use the **Invert Image Profile** in the Peak Detection settings.


8.8 Copy/Paste Regions:

You can copy the blue boundary lines for analyzing same bands on different regions using **Copy Regions** button and then paste using the **Paste Regions** button. Useful for comparing recovery/quantity.

Click **OK** to accept the calculated areas and close the dialog.

8.9 Viewing & Exporting Results

Click the **View Full Results and History** button to open the Analysis Results window.

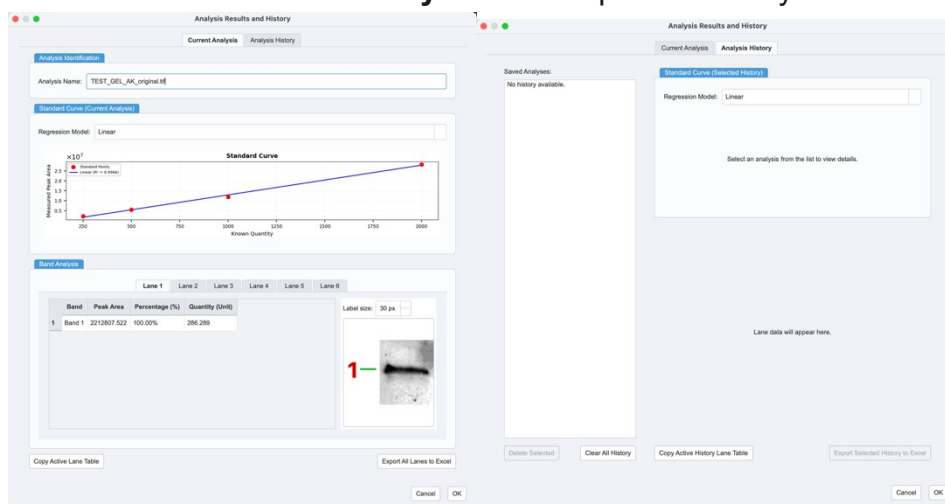


- **Current Analysis** tab: Shows a detailed table and a plot of the standard curve for your most recent analysis. You can copy the table or export all data to an Excel file.
- **Analysis History** tab: Every time you click "OK" in this window, the current analysis is saved to a history file. You can browse past results, view their data, and export them to Excel.

Exporting to PDF: The PDF report bundles an annotated image of your gel (with clearly-labelled lane-number badges), the full parameter and method summary (justified for readability), a data-quality / band-saturation check, per-lane band images, and the results tables.

9. Molecular Weight Prediction

This feature allows you to estimate the size of a protein band based on your standard markers.

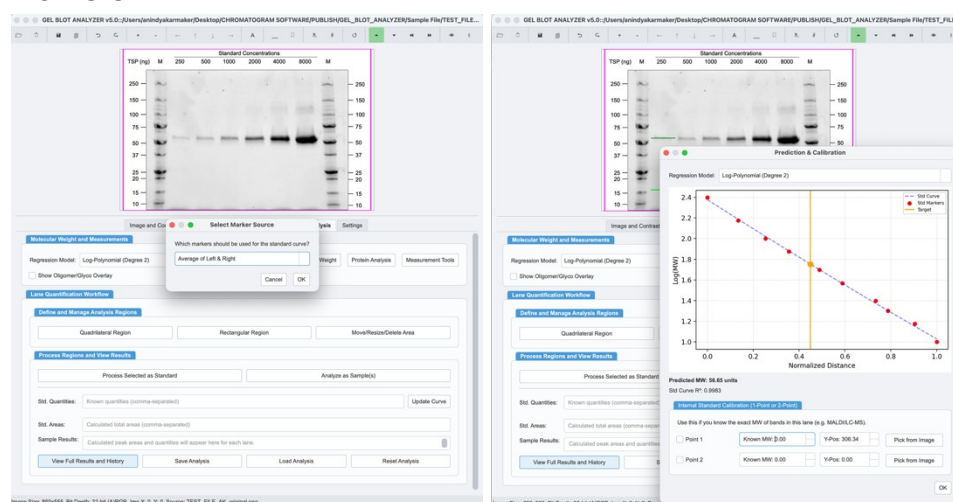

9.1 Setting up the Standard Curve

1. Place your standard markers using the tools in the **Markers** tab. You need at least two numeric markers on either the Left or Right side.
2. Go to the **Analysis** tab.
3. Click the **Predict Molecular Weight** button (**Ctrl/CMD + P**).
4. If you have markers on both sides, you will be asked which set to use (Left, Right, or an Average).
5. The cursor will change. Click on the center of the band you want to analyze.
6. The **Prediction and Calibration** dialog appears.

This window shows a plot of your standard markers (log scale) and the regression line used for the calculation. The predicted molecular weight and the R^2 value (goodness of fit) are displayed. You can change the Regression Model to see which one best fits your standards. The available models are Log-Linear (Degree 1), Log-Polynomial (Degree 2), Log-Polynomial (Degree 3), and Log 4-PL.

NOTE: Usually **Log-Polynomial (Degree 2)** fits best for a gradient gel for accurate estimation without overfitting.

Click **OK** to accept the prediction.

10. Advanced Features

10.1 Protein Size Analysis (Oligomers & Glycosylation)

This tool helps you hypothesize the identity of bands based on a protein's primary sequence and potential modifications.

The screenshot displays the 'Protein Size Analysis' software interface. The main window is titled 'Protein Size Analysis' and contains several panels. The 'Protein Sequence Input' panel on the left has a text box with a protein sequence and an 'Analyze Sequence' button. Below it, the 'Sequence Analysis Results' panel shows 'Total Potential Sites Found: 3' and 'Calculated Potential Fragments' with a list of molecular weights. The 'Physicochemical Properties' panel in the center lists values for Isoelectric Point (pI), Sequence Length, Ext. Coeff. (Reduced), Absorbance (Reduced, 0.1%), Ext. Coeff. (Oxidized), Absorbance (Oxidized, 0.1%), and Net Charge at pH. The 'Fragment Modeling Parameters' panel below that includes fields for Base Protein MW (Da), Avg. Glycan Mass (kDa), Number of Glycans, and Number of Oligomers. On the right, the 'Analysis' tab is active, showing a gel image with molecular weight markers and a 'Molecular Weight and Measurements' section with regression model, quantification workflow, and analysis results. At the bottom, there are buttons for 'Load Analysis', 'Export Analysis', 'Cancel', and 'OK'.

1. Click the **Protein Analysis** button on the **Analysis** tab.
2. Paste your protein's amino acid sequence into the top text box.
3. Click the **Analyze Sequence** button. This will automatically calculate the theoretical **Molecular Weight (MW)** and populate the "Base Protein MW (kDa)" field. Moreover, it will calculate the **Molar Extinction Coefficient** (for both reduced and oxidized forms) based on Tryptophan, Tyrosine, and Cysteine content, identify potential **N-Glycosylation sites** (N-X-S/T sequons) and display a formatted sequence view with residue numbering and a summary of the analysis.
4. In the dialog, you can now fine-tune the model: **Base MW:** The calculated value can be manually overridden if you have a more accurate experimental value. **Avg. Glycan Mass:** Select a common glycan type from the dropdown menu to auto-fill the mass, or select "Custom..." to enter your own value (e.g., 2.5 kDa). **Number of Glycans:** This defaults to the number of potential sites found but can be changed to model different glycosylation states. **Number of Oligomers:** Set the number of oligomeric states to model (e.g., 2 for dimer). The **Calculated Potential Fragments** box will automatically update to list the theoretical molecular weights for all combinations of the parameters you have set.
5. You can save or load your entire analysis using the buttons at the bottom:
 - **Export Analysis:** Saves all the information (sequence, parameters, results) to a readable text file.
 - **Load Analysis:** Loads a previously saved analysis file, repopulating all fields in the dialog.

6. After clicking **OK**, if you have an active MW prediction model, you can check the **Show Oligomer/Glyco Overlay** and click **Predict Molecular Weight** button to show the overlay as lines on your image.

Show Oligomer/Glyco Overlay: When enabled (after a successful Molecular Weight prediction), the software overlays the predicted masses of every modeled oligomer and glycosylation state directly onto your gel. Each predicted band is drawn from a reference position with a colored, pointed arrowhead and a matching-colored label showing its mass in kDa.


The labels are grouped by oligomer: every form of one oligomer (its base mass plus all of its glycosylated variants) shares one aligned column, and each successive oligomer group steps to the right in a clean cascade that begins at the end of the previous group, so labels stay readable even when many bands are close together.

Show Partial Glycosylation: This checkbox controls how glycosylation states are listed. When it is OFF, only the fully-glycosylated form of each oligomer is shown. When it is ON, every partial-occupancy band is shown, including the unglycosylated backbone. For example, a 50 kDa protein with three 2 kDa glycan sites shows 50, 52, 54 and 56 kDa with the option ON, and only 56 kDa with it OFF.

For oligomers, the number of glycosylation sites scales with the number of subunits — a dimer of a protein with 3 sites has 6 sites in total — so both its fully-glycosylated mass and its partial ladder are calculated accordingly.

10.2 Automatic Lane Detection


This feature helps you quickly place markers for a set of straight or skewed lanes.

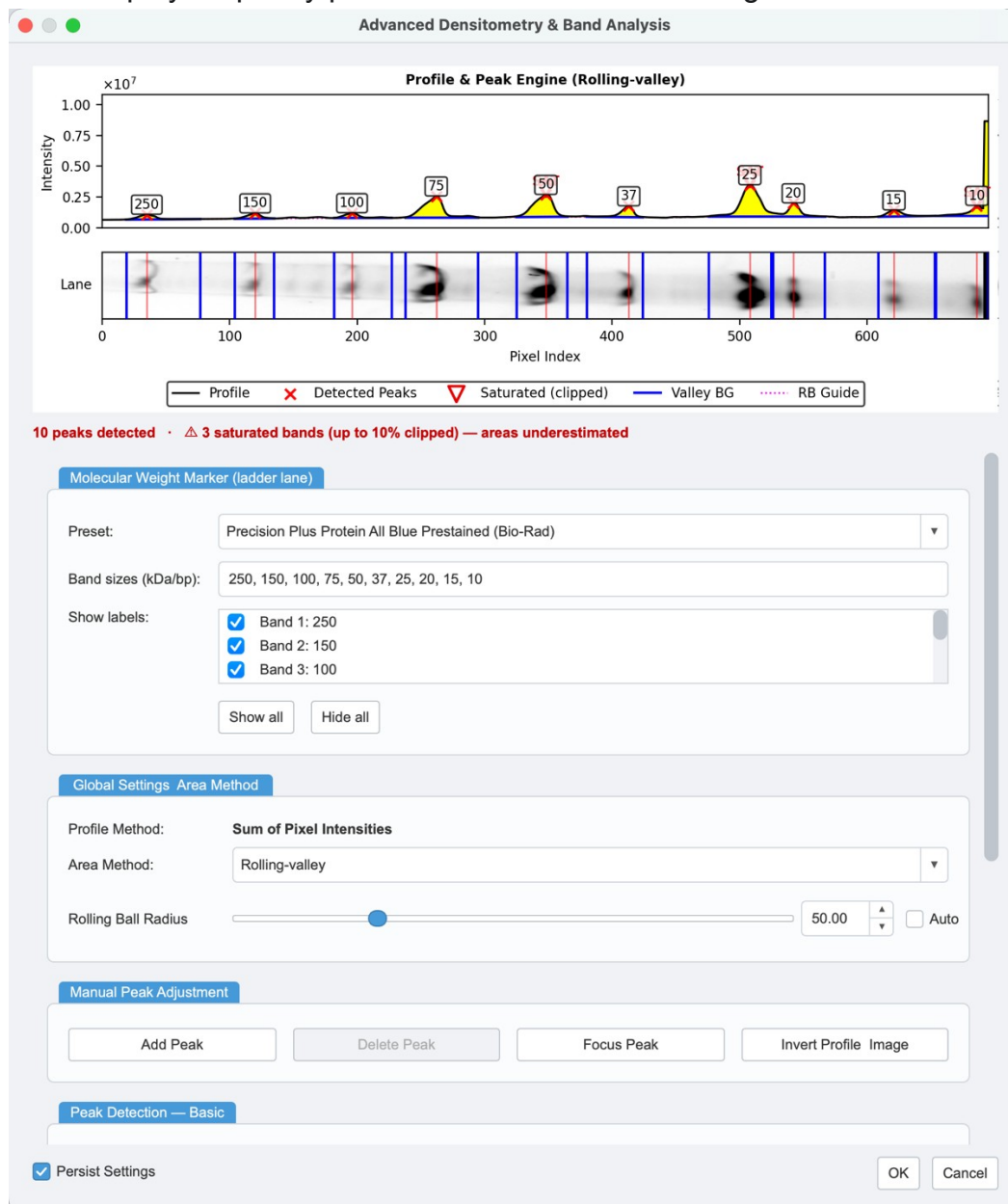

1. Click the **Automatic Lane Markers** button on the toolbar.
2. **Select Marker Side** (Left or Right) to place the markers.
3. **Select Lane Region Type** ("Rectangle" for straight lanes, "Quadrilateral" for skewed).
4. Draw the region covering all the bands in your lane.
5. The **Advanced Densitometry & Band Analysis** dialog will open, showing the intensity profile. Adjust the sliders to ensure all your bands are detected as peaks ('x' markers). Use the delete key to delete unintended peaks

6. You can select the default Preset options or use the **Custom** preset to put the band sizes. You can also hide band sizes by unchecking from **Show labels**.
7. Click **OK**. The software will place the standard markers at the position of each detected peak.

10.3 Automated One-Step Gel Analysis (Auto Analyze Gel)

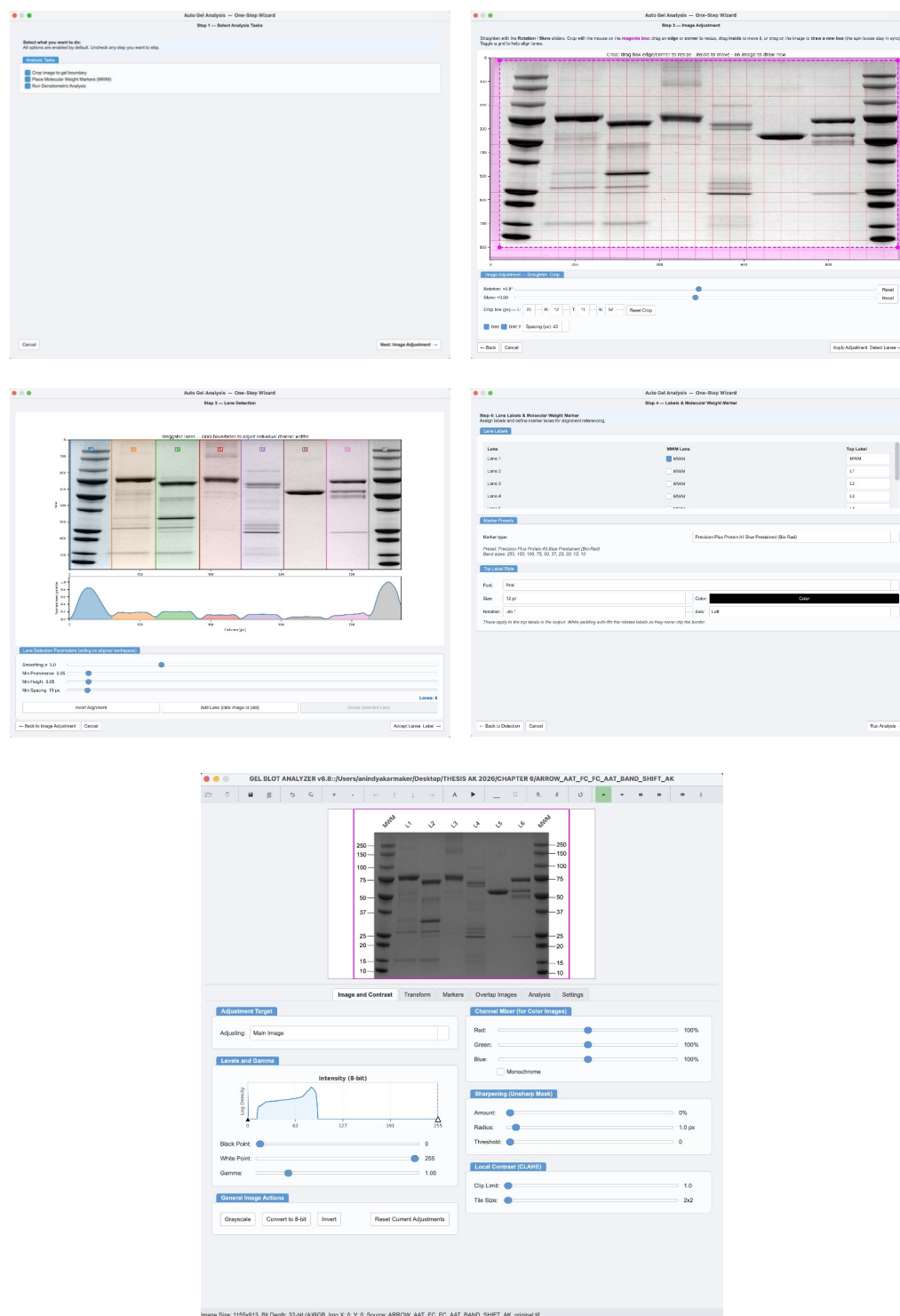

The Auto Analyze Gel tool runs an entire workflow—from a raw image to densitometry results—through a single guided wizard. Launch it from the “Auto Analyze Gel” button on the toolbar or from the Tools menu. It is ideal for quickly quantifying an entire gel without manually defining each lane.

When you start the tool, you are first asked how to draw the gel boundary: choose Rectangle for a straight gel or Quadrilateral for a skewed/tilted gel, then draw the region around all of your lanes. The One-Step Wizard then walks you through five phases:

1. **Task Selection:** Choose whether to run full densitometry with quantification or a marker/labeling-only pass.
2. **Image Adjustment:** Straighten the gel with the Rotation and Skew sliders, fine-tune the crop with a draggable crop box, and toggle a grid for alignment.
3. **Lane Detection:** The software builds a column-intensity profile and automatically places lane boundaries on the prepared canvas. Drag the boundary lines to add, remove, or reposition lanes. **Resize the lanes first** and then add missing lanes. It will number them accurately.
4. **Labels & Molecular Weight lane:** Pick the molecular weight marker lane and rename the lane top-labels
5. **Band Detection & Analysis Type:** Adjust the molecular weight marker options and adjust the band/peak detection settings. For peaks not detected, can manually add them or adjust the height/prominence settings

When you finish, the wizard applies all geometry (rotation, skew, crop, padding) to a full-color, full-resolution copy of your image, places the markers, runs densitometry on every detected lane, and opens the Analysis Results window with the combined data.

10.4 Measurement Tools

The Measurement Tools window—opened with the Measurement Tools button on the Analysis tab—provides image calibration and on-image measurements in a separate, always-on-top window, so you can keep working in the main viewer. It offers three tools:

1. **Set Scale:** Calibrate the image by drawing a line across an object of known length (e.g., a ruler or printed scale bar) and entering its real-world length. The active scale is shown in the window.
2. **Measure Distance:** Click two points to measure the straight-line distance between them. The result is reported in calibrated units once a scale is set, otherwise in pixels.
3. **Measure Area:** Click and drag to draw a freehand shape and calculate its enclosed area.

Use Clear All to remove the current scale and any measurement overlays. Closing the window deactivates the active measurement tool.

10.5 Application Settings (Settings Tab)

The Settings tab (the sixth tab) holds application-wide preferences. Click Apply and Save Settings to store them. Some options take effect immediately, while others require restarting the program.

1. **Processing Unit Priority:** Choose whether image processing runs on the GPU or the CPU. When GPU is selected you can pick the target OpenCL device (GPU:0 is typically the integrated GPU, GPU:1 the discrete GPU). A live label shows the

currently active device. Switching CPU/GPU applies immediately; changing the specific GPU ID requires a restart.

2. **User Interface Scaling:** Adjust the overall UI scale (50%–250%) to suit high-DPI or small screens. Requires a restart.
3. **Application Dimensions:** Set the image-viewer width and height and the program UI width (in pixels); changes take effect after a restart. The viewer now scales its content to these dimensions — the image, markers and labels all zoom together to fill the viewer. Larger values simply give a bigger, proportionally-scaled view.
4. **Reset All Settings:** Restores all application preferences (viewer size, UI scale, CPU/GPU choice and saved marker presets) to their defaults by deleting the saved configuration file. You are asked to confirm first. This does NOT affect your per-image saved configuration files. A restart is recommended afterwards.

11. Keyboard Shortcuts

Action	Windows/Linux	macOS
File		
Load Image	Ctrl+O	Cmd+O
Save Image & Config	Ctrl+S	Cmd+S
Edit		
Undo	Ctrl+Z	Cmd+Z
Redo	Ctrl+Y	Cmd+Shift+Z
Copy Rendered Image	Ctrl+C	Cmd+C
Paste Image	Ctrl+V	Cmd+V
View		
Zoom In	Ctrl+=	Cmd+=
Zoom Out	Ctrl+-	Cmd+-
Pan View	Arrow Keys	Arrow Keys
Tools & Modes		
Reset Image & Annotations	Ctrl+R	Cmd+R
Predict Molecular Weight	Ctrl+P	Cmd+P
Clear Analysis Data	Ctrl+Shift+P	Cmd+Shift+P
Place Left Marker	Ctrl+Shift+L	Cmd+Shift+L
Place Right Marker	Ctrl+Shift+R	Cmd+Shift+R
Place Top Marker	Ctrl+Shift+T	Cmd+Shift+T
Place Arrow Marker	Ctrl + Arrow	Cmd + Arrow
Toggle Grid and Snap (X & Y)	Ctrl+Shift+N	Cmd+Shift+N
Toggle Line X	Ctrl+Shift+X	Cmd+Shift+X

Toggle Line Y	Ctrl+Shift+Y	Cmd+Shift+Y
Increase Grid Size	Ctrl+Shift+ ↑	Cmd+Shift+ ↑
Decrease Grid Size	Ctrl+Shift+ ↓	Cmd+Shift+ ↓
Toggle Guide Lines	Ctrl+G	Cmd+G
Convert to Grayscale	Ctrl+B	Cmd+B
Invert Image	Ctrl+I	Cmd+I
Toggle Snap to Grid	Ctrl+Shift+N	Cmd+Shift+N
Toggle Dark/Light Mode	Ctrl+Shift+D	Cmd+Shift+D
Quickly Switch Tab	Ctrl+1/2/3/4/5/6,	Cmd+1/2/3/4/5/6,
Exit Any Active Mode	Esc	Esc

12. Troubleshooting

- **"A required library is missing" error on startup:** Make sure you have installed all the dependencies listed in the [Installation](#) section.
- **Molecular Weight Prediction is grayed out:** You need to place at least two standard markers on the Left or Right side before this option becomes available.
- **Image looks distorted after rotation or skew:** These operations are previews. Click the "Apply" button in the **Transform** tab to make the change permanent and see a high-quality result.
- **Cannot move or resize a lane area:** Make sure you have clicked the **Move/Resize Selected Area** button on the **Analysis** tab to enter the correct mode.
- **The application crashes unexpectedly:** An error_log.txt file is created in the same folder as the application (next to the executable). It records the full error details for every kind of failure — including errors raised inside the interface, in background tasks, and low-level (native) crashes — so it is the first place to look when reporting a problem. Each unexpected-error dialog also has an expandable Details section containing the same traceback.